

ASEPTIQUIK® W ASSEMBLY AND BEST PRACTICES

ASSEMBLY STEPS

1 FLIP

FLIP the pull tab down while disengaging it from the shipping clip until fully open and clicked into place.



BEST PRACTICES




Tips: Press down on the shipping clip latch as you flip the pull tab.

Note: When the pull tab is opened, the bar weld will peel away as intended. The sterile boundary created by the circular weld remains intact throughout the process.


2 CLICK

Align the two halves and press together until you hear a **“CLICK”** sound of the 1st stage connection latches. Then, Click the pull tabs together.

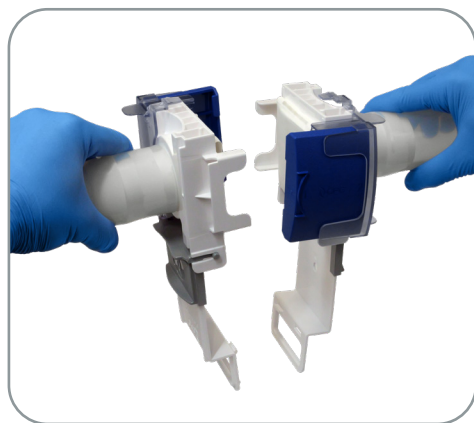


Tip: Some maneuvering may be required to slip the body hooks into the body latches. Multiple hand placement options – see images on the right.

Hand Placement Option 1

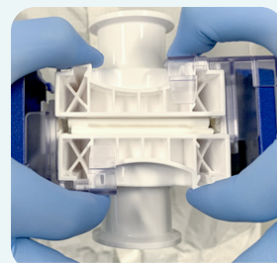

Hand Placement Option 2

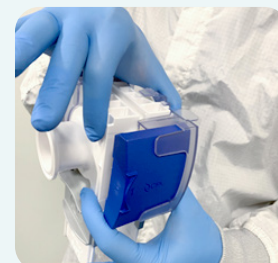

Hand Placement Option 3



Correct Engagement



Incorrect Engagement



Final Step: Click tabs



3 PULL

PULL the sterile barrier components by placing one hand in the handle of the pull tabs and the other hand at the hose barb or sanitary flange termination.



Side Latch (blue) Shipping Clip (clear)



Tips: Completely remove pull tabs and sterile barriers by pulling directly away from connector body. Avoid squeezing the pull tab guides inwards during the sterile barrier pull process.

IMPORTANT: Do not proceed to step 4 until the pull tabs and sterile barriers are fully removed from the connector.

Pull Tab Guide (Gray)

4 SNAP

Remove the clear shipping clips, then **SNAP** blue side latches “1” first, followed by side “2” latches, and confirm they are attached to the body.



Remove shipping clips.



IMPORTANT: Side 1 latch must be clicked prior to Side 2 latch, to ensure proper connection.

